


HOMEPRO
JFK Health

BADY MASSAGER

JFK-896B

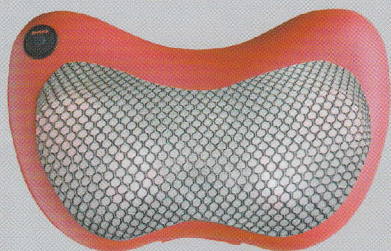

HOMEPRO

JFK Health

JFK-896B

INTRODUCTION

The Intelligent Multi-function Kneading Massager use the high-tech computer digital technology, the massage mechanism as massage staff, Integration of the theory of human engineering, promote blood circulation, eliminate the fatigue, the pain of the neck and relieve the working pressure.

CHARACTERISTIC

The product is exquisite and portable, it's not only can be enjoy at home, but also in car and office, give you an unexpected feeling relaxed, so that you work to achieve a multiplier effect.



FUNCTION

- The four heating and kneading massage head, be able to more in-depth massage of the fatigue muscle tissue, promote blood circulation. Alleviate pain, relax. Allows you to more comprehensive and more comfortable to enjoy the pleasure of massage.
- The smooth heat treatment, have hot compression effect.
- Automatic dual-rotating massage function, enjoy two different style massage.
- Automatic time setting for 20 minutes.



Relieves the tightest knots



Stretches and relaxes muscles

OPERATION INSTRUCTION



- **Home use:** Insert adapter's output plug onto the jack on the product, and then plug the adapter onto the electrical outlet box.
- **Car use:** Insert adapter's output plug onto the jack on the product, and then plug the cigar lighter onto the cigar lighter receptacle on car. Get the power cord connected, the product is in standby mode.
- Press the button at the top left corner, the product start to work, massage with rotating, kneading and heating, and the kneading massage will automatic change the direction for every minutes.
- Press the button again, stop heating, continue kneading.
- The third press the button to restart the heating, the function of the button circle like this.
- Press the button for 2 seconds, the product stop to work.
- To prolong the durability of the massager, it have automatic time setting for 20 minutes, wait half an hour before reusing.



ATTENTION

- An adapter with wall outlet box and car outlet is provided in each unit, do not use any other facilities for power source.
- stop using the product whenever any defects was founded, or the product works unsmooth. Send it to the dealers and do not open it.
- The product should not be used for more than 20 mins.
- Do not draw or pull the power cord after use, but gently hold the end of the power cord and unplug.
- Stop using the product whenever any defects were found on the skin of power cord, contact your dealer to change it.
- Do not use the product in a wet or moist condition.
- The product should not be used by or near children.

CAUTION

If you are under medical treatment or are included in any of the following groups below, consult your physician and follow his or her directions regarding usage of this massager: heart disease, malignant tumor, pregnancy and menstruation, the back deformity, have the history of treatment use in the affected part, osteoporosis, implant heart pacemaker, forbident to use.

ELECTRIC PARAMETER:

- Name: Body Massager
- Modle: JFK-896B
- Rated voltage: AC100-240V / DC12V
- Rated frequency: 50-60Hz
- Rated power: 24W
- Rated using time: 20Mins
- Safety standard:



BADY MASSAGER
JFK-896B