



Thank you for using our product. Please read the user guide before use and keep it for future reference.

Introduction	2
Features of the product	2-3
Preparations before use	4
Names and Functions of Parts	5
Use Methods	5-7
Security Matters	7
Overhaul and Specifications	8
Product Maintenance	9


Instruction

Product principle

Based on techniques of far-infrared heating and micro computer chip control, the eye massager brings the following effects on the eyes by heating, vibrating and massaging the eye area: it increased the oxygen content around the eye area, improves the adjusting ability of the crystalline lens, recovers the elasticity of ciliary muscles, improves the blood circulation around the eye area, relieves visual stress, reduces eye fatigue, reduces visual disturbance, stimulates metabolism and maintains visual health.

Application scope

The product can function as an assistant treatment for pseudomyopia, mild myopia and amblyopia, and take care of health by reducing eye fatigue and preventing myopia.

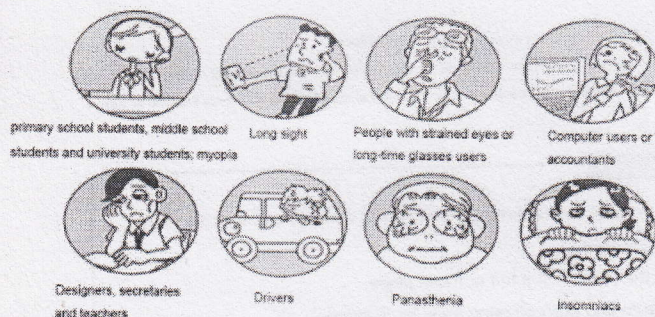

Features of the product

The product is a mixed mode high-tech eye massager.


The product uses humanized analog control during the whole process by adopting computer chip control technology.

The product adopts world's leading high-tech achievement: infrared heater. The material is characterized by even heating, low power consumption, high density and intensity of far infrared rays. Infrared rays can activate cells exposed to radiation, cause molecular resonance, strengthen the binding force between molecules, improve the regeneration capacity of tissues, enhance immunity, increase oxygen content and stimulate blood circulation, through which health is taken care of.

The product successfully applies the brand new concept of air pressure to eye care, which humanizes the product. Under the control of the micro computer, air pressure conducts all-around humanized massaging in the eye area, relieves eye pressure quickly, reduces eye fatigue, recovers the contractility of ciliary muscles, adjusts the biological elasticity of ciliary muscles and prevents the ocular axial from elongation or atrophy.

The product adopts embedded three-dimensional magnetic field; during use, the magnetic field in constant motion acts upon the main facial acupuncture points and capillary vessels around the eyes by cutting magnetic lines of force, thereby generating micro current inside the capillary vessels which clears and activates the channels and collaterals and harmonizes qi and blood.

Functions of acupuncture points

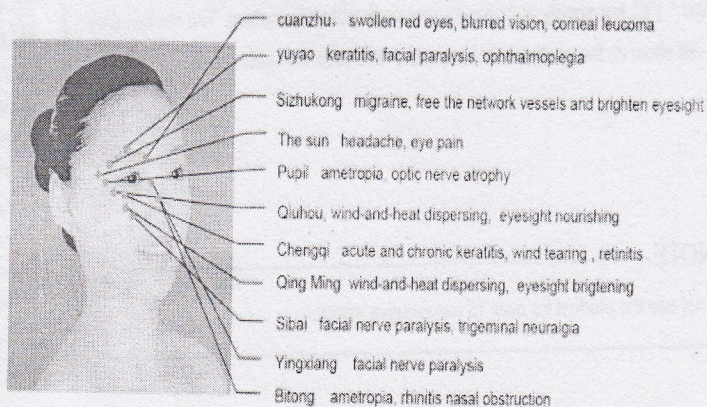




Preparations before use

Installation of batteries

1. Press the battery cover slightly by hands and push it outwards to release it.
2. Take out the hand controller battery cover and fit in 4 No.5 batteries.
3. Close the battery cover.

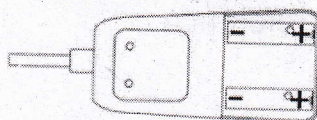


NOTE

- Confirm that the anodes and cathodes of the batteries are in right place and the battery case cover is full closed.
- As for the anode and cathode of the battery case, please refer to the picture above; alkaline batteries are recommended.

Operation methods of AC adaptor

1. Please use DC6V stabilized voltage supply allocated by our company.
2. Connect the AC adaptor to the power input socket of the hand controller.
3. Connect the AC adaptor to the power supply socket.

How to put on the product

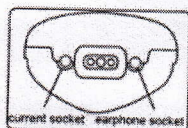

1. Earphone connection
Connect the earphone to the earphone socket of the product.



2. Put on the product
Please put on the product with the cloth cover in position of the eyes, and adjust the belt and bucket it for the most



3. Adjust the cloth cover for the most comfortable position.

14/

Self-selected mode

Press the "□" key and choose your preferable massage mode.

Vibration massage mode

When using the vibration massage mode, press the "☺" key and shift between ten different vibration massage modes.

MP3 music procedure

Press the "▶/||" to play or pause music, and press the "◀◀" key or the "▶▶" key to change music.

A USB data line can connect the product to the PC computer for downloading and charging.

Time setting

The product is provided with three time settings: 5 minutes, 10 minutes and 15 minutes.


Press the "⏸" key shortly and choose your preferable time setting. The chosen time setting will show on the controller (the intended time area will lighten up).

NOTE

- Do not use the product for over 15 consecutive minutes.



Names and Functions of Parts

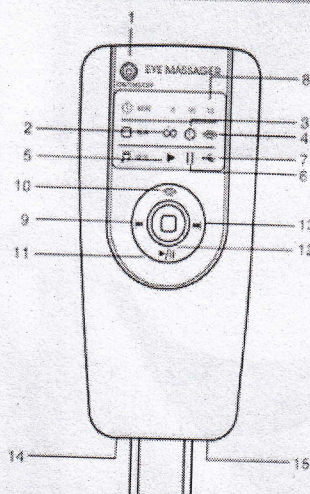

1. On/off/time key (a short press for turning on the product, another press to set time, and a long press for turning off the product)
2. Pressure mode indicator light
3. Hot pack mode indicator light
4. Vibration mode indicator light
5. MP3 play indicator light
6. MP3 pause indicator light
7. USB computer downloading indicator light
8. Time setting indicator light
9. Previous song
10. Vibration mode selection key
11. MP3 play/pause
12. Massage mode selection key
13. Next song
14. DC interface (regulator specifications: DC6V 10000MA)
15. Earphone interface



Use Methods

Press "⏸" the button to start up

The product will automatically goes into the preset massage procedures, and the pressure indicator light "∞" will lighten up. If the user needs to stop the ongoing massaging, press the "⏸" button and hold the pressing.

Massage procedures (one-key circulation design)

Press the "□" button and choose from three modes.

Relaxation procedure: the pressure mode indicator light "∞" of the controller will be on.

Lively procedure: the pressure mode indicator light "∞" and hot pack mode indicator light "☺" of the controller will be on.

Quality mode: the pressure mode indicator light "∞", the hot pack mode indicator light "☺" and the vibration mode indicator light "☺" of the controller will be on.

After use

After use, press the "⏸" key to turn off the controller; if the product is not to be used, turn off the product, unplug the AC adaptor from the socket and take out the batteries.



Security Matters

1. Please read carefully and understand the whole content of the user guide, and operate the product accordingly.
2. People with symptoms of eye disease are not allowed to use the product; patients are forbidden to use the product who had eye operations or have serious eye diseases such as glaucoma, cataract, retinal detachment, etc.
3. Before use, take off the contact lens, hairpins, hair ornaments, etc.
4. It is forbidden to clean the product with water; please use alcohol-contained cotton to clean the product.
5. It's forbidden to disassemble the product without authorization.
6. The attachment of the air tube on the right side of the product makes the weight of the product on the head feel imbalanced, which is normal and won't cause negative effect on the use of the product.

NOTES

1. Put on the product before turning it on to avoid damage of the air sac from excessive pressure.
2. Do not use the product in a hot or moist place.
3. Strong sunlight and gas corrosion must be avoided.
4. Do not use the product during outdoor walking or driving.
5. Handle the product gently and do not soak it in liquid.
6. Please use No.5 batteries or the allocated voltage regulator. When exterior power supply is in use, it is forbidden to put batteries inside the battery case. When the product is not used, please take out the batteries to avoid malfunction caused by battery liquid leakage.
7. Please use the product with eyes closed. During use, the vision becomes vague at first and gets clearer gradually after closing the eyes for 1-3 minute.

Overhaul and Specifications

If any problem occurs during product operation, please go through the following examination. If the problem still can not be solved or the product still can not work properly, please turn off the product and make contact with our service center for product repair.

Please do not disassemble or repair the product without authorization.

Problem	Possible cause	Treatment
Failure to start	Is the battery power deficient?	If yes, please change the batteries.
	Is the battery power deficient?	Make sure the anodes and cathodes of the batteries are in right position.
		Please fix the battery in the battery box
	Is the voltage regulator proper connected?	Make sure the connection is right.
	Are the electric wires intact?	Examine the connection of the wires; contact the manufacturer if any damage occurs.
The indicator light is dim. The air pressure is too weak.	Is the batter power too low?	If yes, please change the batteries.
	Is the batter power too low?	If yes, please change the batteries.
		Make sure the batteries are installed and fixed in the battery case.

Technical Specifications	
Work environment	Temperature 5℃ -40℃ ,relative humidity≤ 80%
Atmospheric pressure	860hpa-1060hpa
Magnetic field intensity	50mt-250mt
Hot pack temperature	30℃ -42℃
Power supply	6V DC
Input power	6W
Net weight	485g

For keep

Improper treatment of discards will cause environmental pollution



Product Maintenance

Power supply examination

Examine the power supply and do not turn on the product if the voltage exceeds the specified limit.

Maintenance

If the product is dirty, please use dry soft cloth to wipe it.

If the product is very dirty, use alcohol-contained soft cloth to wipe the body of the product.

Put the product in a safe place after use and do not put weight on the product.

Take ant-moisture, anti-mouse and anti-pollution precautions.

When failure occurs, send for professional technicians for repair and do not disassemble the product without authorization.

Our company reserves rights of improving the product and the user guide and apologizes for not giving extra notice if some change happens. Users should consult the practical product for its looks.