

VitaDental


Electronic toothbrush for a healthy mouth

Operating instructions VitaDental electronic toothbrush – *Your electric toothbrush*

Charge: 2 AAA NiMH batteries

Output: 1 W

Performance: 8000 oscillations per
minute

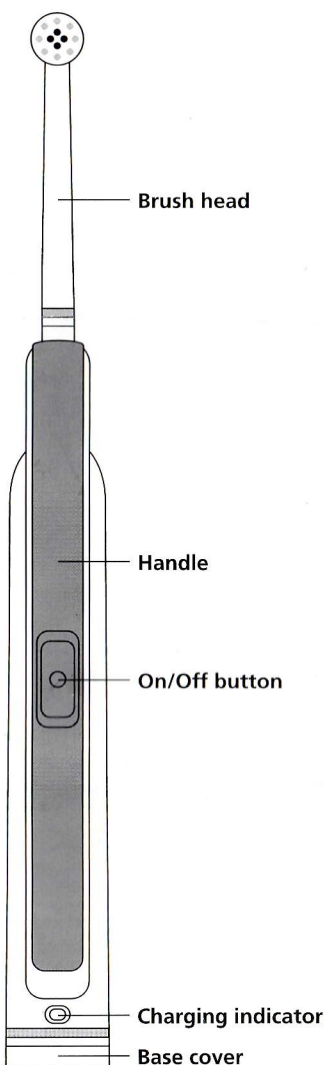

Important

Check your electric toothbrush regularly. Do not use the product if the cord or appliance has been damaged.

- ✓ The product should not be used by children under the age of 3.
 - ✓ Do not place or store the charger so it may be at risk from falling or getting pulled into a bathtub or hand basin.
 - ✓ Do not place the charger in water or any other liquid.
 - ✓ Do not try to pick up the charger if it has fallen into water. Unplug it first.
 - ✓ Do not use AAA alkaline batteries with the toothbrush if you intend to charge it. This may cause an explosion.
 - ✓ Product is only for intended use and in accordance with these operating instructions.
 - ✓ The batteries are recyclable.
-

Features

- ✓ Oscillates clockwise and counter clockwise 8000 times per minute.
 - ✓ Remains charged for up to 14 days when brushing your teeth for two minutes, twice a day.
 - ✓ To get the best results and optimal dental hygiene the brush head should be replaced every month.
-

Using the electric toothbrush

You can use any toothpaste.

Dip the brush head in water and add toothpaste. Avoid splashing by placing the brush head against your teeth before turning the toothbrush on. Brush your teeth for approx. 2 minutes. Move the brush from tooth to tooth and hold still on each tooth surface for a couple of seconds. Brush both gums and teeth.

Clean the brush head when you are finished. Change the brush head by pulling it straight up and replacing it with a new one.

Maintenance

- ✓ Make sure you rinse the brush head properly after use.
- ✓ Your electric toothbrush contains rechargeable NiMH type batteries. Do not charge the batteries until they are completely discharged, otherwise the battery cells may become damaged.